




# VirusID Pro

## AI-Driven, Rapid, Portable Multi-Virus Detection Kit

# Technology & Process

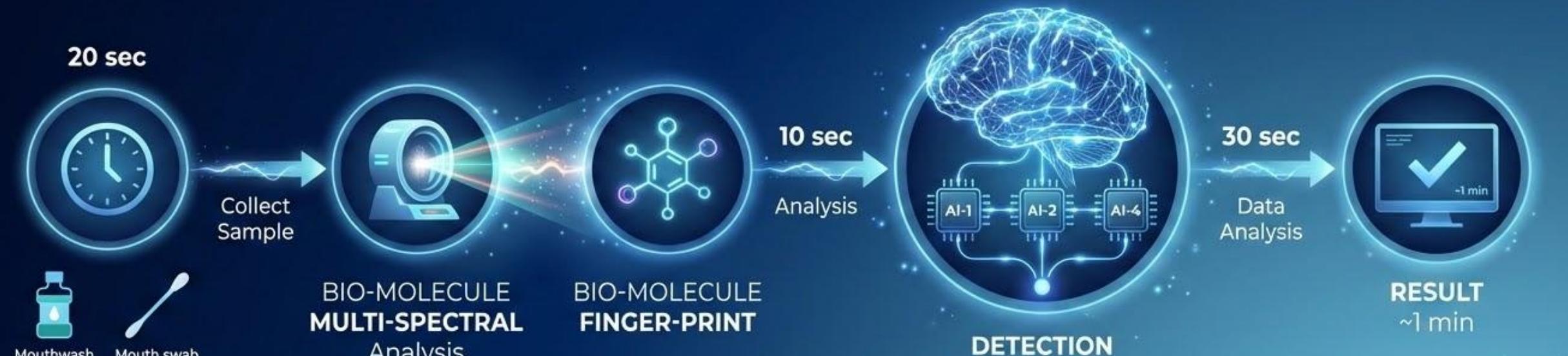




Current Technologies vs VirusID Pro	
• Indirect diagnosis	✓ Direct Multi-Virus diagnosis from one sample
• Test Kits, consumables, logistics issues	✓ No Test Kits
• Multiple Equipment	✓ Only one test Equipment
• Laboratories, medical staff and facility (for PCR based tests)	✓ No medical staff
• Time (more than 15 minutes) & Cost	✓ Time (~1 minute) & Lower cost

## PoC Study 1

Can VirusID Pro be trained to detect and differentiate between various biomolecules, such as Denosumab, the Hepatitis-A vaccine, and the Influenza vaccine, when they are mixed in saliva?

### Test Samples

- Hepatitis A Vaccine
- Denosumab
- Influenza Vaccine

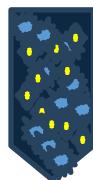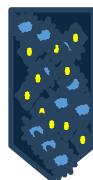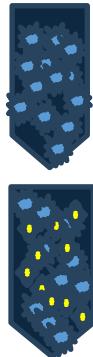

**30 sec**

Measure in VirusID Pro

### RESULTS



+



-



-



+



+



-



+



+



+

This study showed that VirusID Pro effectively learned to detect and distinguish biomolecules and viruses under laboratory conditions

## PoC Study 2

Can VirusID Pro and point-of-care device detect the SARS-CoV-2 virus in mouthwash samples with acceptable accuracy?


### Test Samples

mouthwash samples



mouthwash samples

+  
SARS-CoV-2

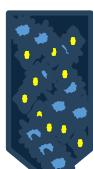

**30 sec**

Measure in VirusID Pro

### VirusID Pro results



-



+

### PCR results



-




+

This study demonstrated the performance of the AI's learning capabilities in heterogeneous mouth saliva samples within a pre-clinical study setting.

# Technology

- **Mobile hardware** collects convoluted/overlapped data obtained via detectors (c1., c2., c3) which are customized to measure the electromagnetic spectrum produced from the interaction of the signals sent from electromagnetic sources (b1., b2., b3 ) to the sample (a). and sends it to the AI-based software.
- In A.I.-based software, a workflow is developed by sequential use of different methods. Then, the data is processed by a few sequential learners.
- The formed raw output data goes to a decision-making algorithm to connect scientific method development to A.I. and finally, the output is given to the user.

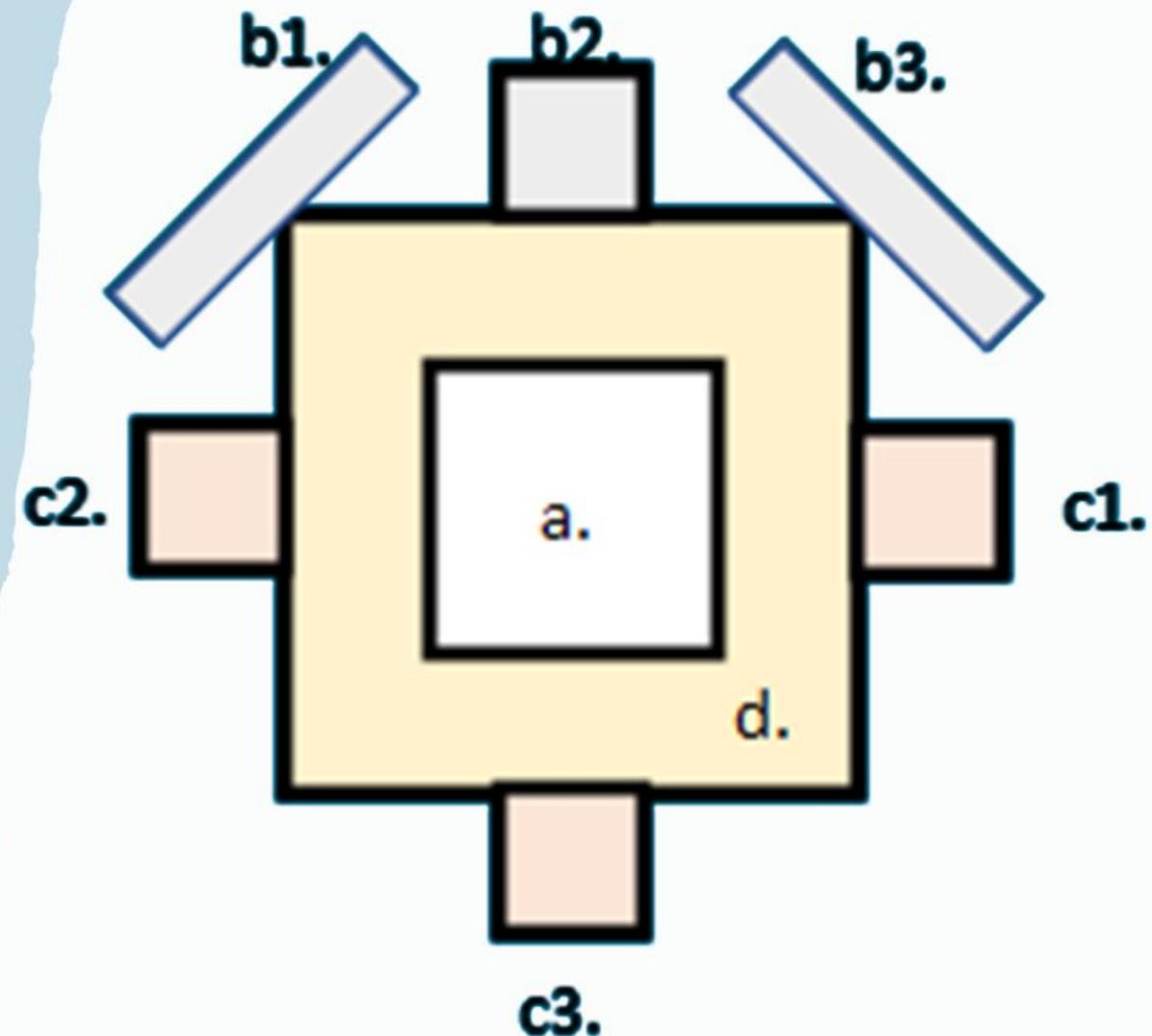

## Clinical Study 1

Is it feasible to develop VirusID Pro for precise COVID-19 detection using nasal swabs and RT-PCR validation?

**RESULTS:** Our findings indicate that VirusID Pro offers a reliable and highly sensitive solution for detecting COVID-19. Its ease of use, rapid results, and suitability for point-of-care testing position it as a promising tool for any infectious agent diagnosis

### BEFORE ATTENTION ADDED INTO MODEL



**2448**

patient samples for training




**816**

patient samples for evaluation




Open Learning study

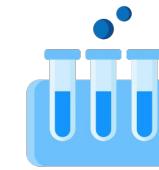

**350**

patient samples

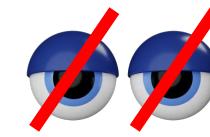

**Double-blind**  
testing



Simultaneous real-life case study

Table-1: Sensitivity, specificity, PPV, and NPV values obtained from the AI learning validation stage with respect to RT-PCR

	PCR (+)	PCR (-)	Total	Sensitivity (%)	Specificity (%)	NPV (%)	PPV (%)
MSA (+)	328	50	378	88.4	88.76	90.18	86.77
MSA (-)	43	395	445				
Total	371	445	816				

MSA: multi-spectrum analysis, RT-PCR: reverse transcription-polymerase chain reaction, NPV: negative predictive value, PPV: positive predictive value

## Clinical Study 2

Is our technology as sensitive as RT-PCR at Sancaktepe Hospital for COVID detection?

Table-3: Sensitivity, specificity, PPV, and NPV values obtained from the AI testing stage with respect to PCR

	PCR (+)	PCR (-)	Total	Sensitivity (%)	Specificity (%)	NPV (%)	PPV (%)
MSA (+)	115	24	139	81.73	81.99	87.81	75.16
MSA (-)	38	173	211				
Total	153	197	350				

MSA: multi-spectrum analysis, PCR: polymerase chain reaction, NPV: negative predictive value, PPV: positive predictive value

### AFTER ATTENTION ADDED INTO MODEL

F1 Score Val: 0.8742


Classification Report Val:

	precision	recall	f1-score	support
0.0	0.88	0.87	0.87	330
1.0	0.87	0.88	0.87	322
accuracy			0.87	652
macro avg	0.87	0.87	0.87	652
weighted avg	0.87	0.87	0.87	652

F1 Score Test: 0.9117

Classification Report Test:

	precision	recall	f1-score	support
0.0	0.94	0.91	0.93	211
1.0	0.88	0.91	0.89	139
accuracy			0.91	350
macro avg	0.91	0.91	0.91	350
weighted avg	0.91	0.91	0.91	350

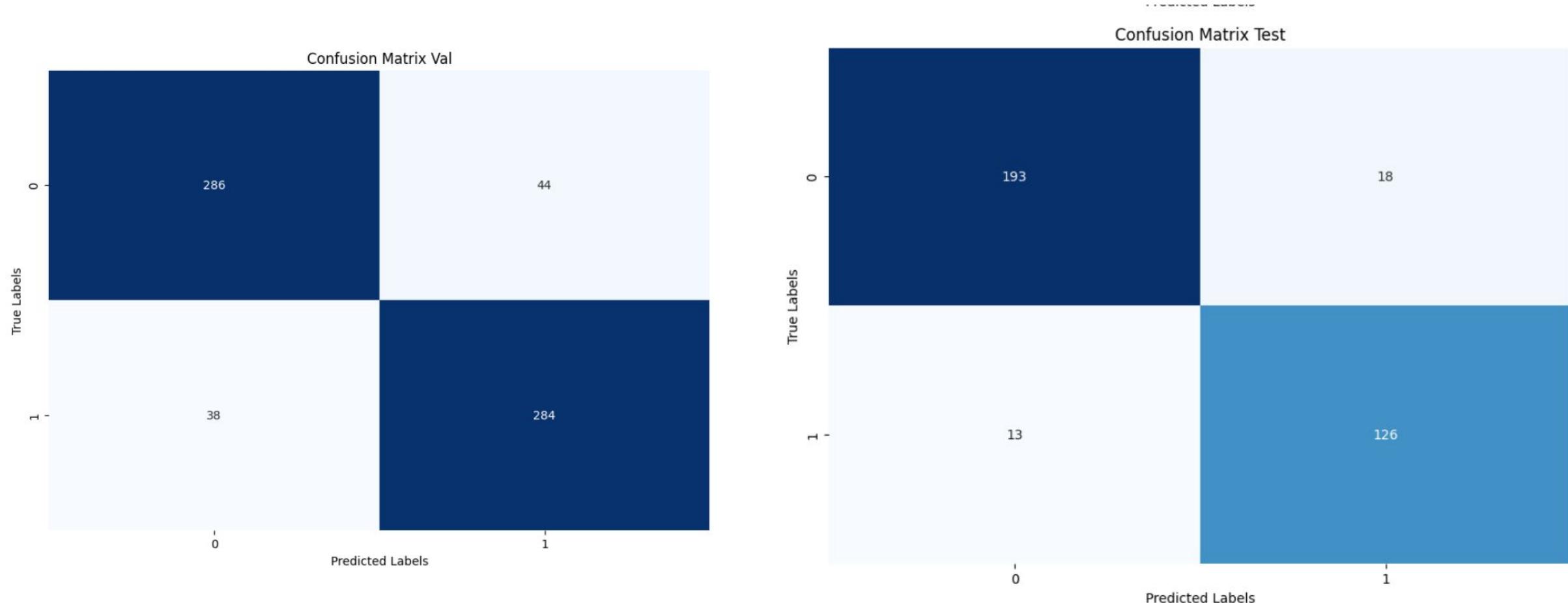

# VirusID Pro – Global Network Vision

Train anywhere & Test anywhere

